GO-TIM

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1. PRODUCT DESCRIPTION AND APPLICATION

Electric wheelchair is designed for people who have problems with mobility, helping disabled patients with their everyday tasks. It can replace a standard wheelchair.

Electric wheelchair is suitable for patients with:

- disability, who have problems with unassisted movement,
- dysfunctions and reduced movement,
- heart and cardiovascular system diseases,
- weakening body due to advanced age.

2. SECURITY INFORMATION

- Read the manual before first use.
- Do not use the wheelchair until fully assembled.
- Use of the wheelchair by people who are unable to operate it independently is not recommended.
- Maximum load capacity must not exceed 136 kg.
- The wheelchair can only be used by one person at a time.
- It is recommended to keep your feet on the footrest and your forearms on the armrests when driving.
- Always check the battery level on the control panel before use.
- This electric wheelchair should only be used on a flat surface or at low elevations. It should not be ridden on a surface inclined at an angle of more than 6 degrees.
- Do not drive the wheelchair at maximum speed when going uphill or downhill.
- Do not use the wheelchair on soft and slippery surfaces such as mud, sand etc.
- Do not ride the wheelchair sideways on highly elevated areas.
- Do not put any additional load such as a trailer etc. on the wheelchair.
- Do not repair or modify the construction of the wheelchair by yourself. Use only original spare parts provided by the manufacturer.
- Do not stand up or sit down on the wheelchair when the control module is switched on.
- Do not use the wheelchair without anti-tip wheels.
- When the freewheel (manual) mode is on, do not get on or off the wheelchair without assistance.
- Do not lean forward to grab items laying in front, on the side or behind the wheelchair. Do not stand on the footrest. It might result in the wheelchair falling over.
- Avoid rapid changes of the direction while driving on elevations.
- Always slow down while driving downhill or making turns.
- Before use check if the wheels are correctly mounted.

- Do not make sudden moves while using the joystick.
- When storing, protect the wheelchair from rain and adverse weather conditions.
- Before turning on the wheelchair, make sure that the levers are set to electric mode.
- Avoid areas with horizontal gaps.
- If the wheelchair is not going to be used for a long period of time – remove the battery.
- Do not use the wheelchair on the road. It might result in a car accident.
- Do not store the wheelchair in places with high temperatures and humidity.

MEDICAL DEVICE INCIDENTS

Any serious incident involving a medical device must be reported to the manufacturer and to the competent authority of the European Union member state where the user or patient resides.

3. GENERAL INFORMATION

Do not operate wheelchair without first reading manual instruction and understanding it. Manual instruction provides user with necessary information regarding safety regulations, maintenance and operation of the vehicle.

Do NOT modify, add, remove or dismantle any of the wheelchair parts or disable functions by yourself. Contact authorized importer service for repairs.

PRE-RIDE SAFETY CHECK

Every time before using the wheelchair, perform inspection of the wearing parts such as brakes, wheels etc.

WEIGHT LIMITATIONS

Do NOT exceed weight capacity (136 kg) of wheelchair as it may result in damage that are not covered by the warranty.

Do NOT carry other passengers on the wheelchair, it may result in instability of a wheelchair or/and accident.

DRIVING ON INCLINES

Do NOT drive wheelchair across an incline or diagonally up or down an incline.

Do NOT drive up or down a potentially hazardous incline (e.g., areas covered with snow, ice, cut grass, or wet leaves).

DRIVING ON DIFFERENT SURFACES

Wheelchair construction allows driving on both outside and inside type of surfaces (concrete, asphalt, floor etc.). However, it is important to remember about few rules to avoid unnecessary risk of damaging the wheelchair.

- Reduce speed while driving on uneven surfaces,

- Avoid driving on grass, mud or loosely packed gravel or sand,
- Do NOT drive on slippery surfaces such as ice as it may result in accident,
- Do NOT drive on surfaces that are unstable,
- Avoid holes and ditches,
- Do not drive on public roads as it may result in vehicle accident. Drive only on sidewalks or special lanes for bikers.

STAIRS AND ESCALATORS

Wheelchair is not designed to travel up or down the stairs or escalator. It can be transported via elevator.

PRECAUTIONS REGARDING WEATHER CONDITIONS

Avoid exposing the wheelchair to inclement weather conditions such as extreme heat or cold, rain, or snow. If you are caught in such conditions while riding, seek shelter and allow the wheelchair to dry completely before continuing or before storing or charging it.

Surface temperature can increase when exposed to external sources of heat (e.g. sunlight).

BATTERY PROTECTION

Always protect the battery from freezing. Do NOT charge battery that has been frozen. If battery is damaged or cracked, immediately contact the importer service and do not use the battery.

Do NOT change any circuit connections. **Any self-repairs are prohibited and might result damages that are not covered by the warranty.**

SYMBOL	MEANING
	Safety Classification type B
	Warning! Indicates a potentially hazardous condition/situation.
	Prohibited! Performing prohibited action can cause personal injury and/or equipment damage. Damage caused by prohibited actions is not covered by warranty.
	Mandatory! These actions must be performed before usage. If user fails to perform the mandatory actions, it may result in personal injury or/and equipment damage

IP54	Protection against limited dust ingress and water spray from any direction.
	Read the manual instruction before operating
	Watch carefully for pedestrians
	Do NOT remove any parts or warning labels of the wheelchair
	Do NOT put fingers or limbs under covers or between movable parts. It may result in injuries
	Do not overload wheelchair. The maximum capacity of the scooter cannot exceed 136 kg
	Do NOT use wheelchair on stairs
	Do not use wheelchair in extreme weather conditions such as heavy rain, extreme humidity, icing etc.
	Warning! Wheelchair surface temperatures may be high
	Use only authorized importer service for repair in case of accident or/and damage to the parts
	Battery contains caustic substances. Dispose according to the local environment regulations
	Before driving, ensure that backrest is in upright position
	Use only lithium maintenance-free sealed batteries from manufacturer

	Keep tools and other metal objects away from the battery terminals. Contact with metal objects may cause an electrical shock
	Disposal and recycling For more information check “Storage and transportation condition” section
	
	The battery charger is for indoor use only
	Flammable materials Do not expose to the open flame
	MR Unsafe
	Consult professionals before using wheelchair in case of mental or physical impairments
	Maintain proper position while driving
	Maintain your focus while driving and avoid distractions
	Fragile
	Protect from humidity and moisture

4. PRODUCT INFORMATION

A. ELEMENTS OF THE WHEELCHAIR

- 1. Handle
- 2. Armrest
- 3. Backrest
- 4. Seat with safety belt
- 5. Control module adjustment knob
- 6. Armrest flip-up button

- 7. Control module with joystick
- 8. Battery & engine connection
- 9. Footrest with calf strap
- 10. Front, swivel 7" wheel
- 11. Rear, solid 11" wheel

- 12. Folding clamp
- 13. Pocket for goods
- 14. Engines and operation mode levers
- 15. Anti-tip wheels
- 16. Lithium battery
- 17. Anti-tip wheel lock

Wheelchair after folding

B. PRODUCT CHARACTERISTICS

- Lightweight & reinforced main frame made of carbon fibre
- Solid rear 11" and front 7" wheels
- Brushless motor
- Control module with LCD display
- Safe and reliable electromagnetic brakes
- Fully closed, maintenance-free battery with handle
- Flip-up armrests, footrests and anti-tip wheels
- Pocket for small goods located at the back
- Two modes of operation: electric and manual

C. TECHNICAL PARAMETERS

Wheelchair cannot be used in areas with flammable gas mixed with air, oxygen or nitrous oxide.

Note! Some parameters might vary depending on user weight, road conditions and battery usage.

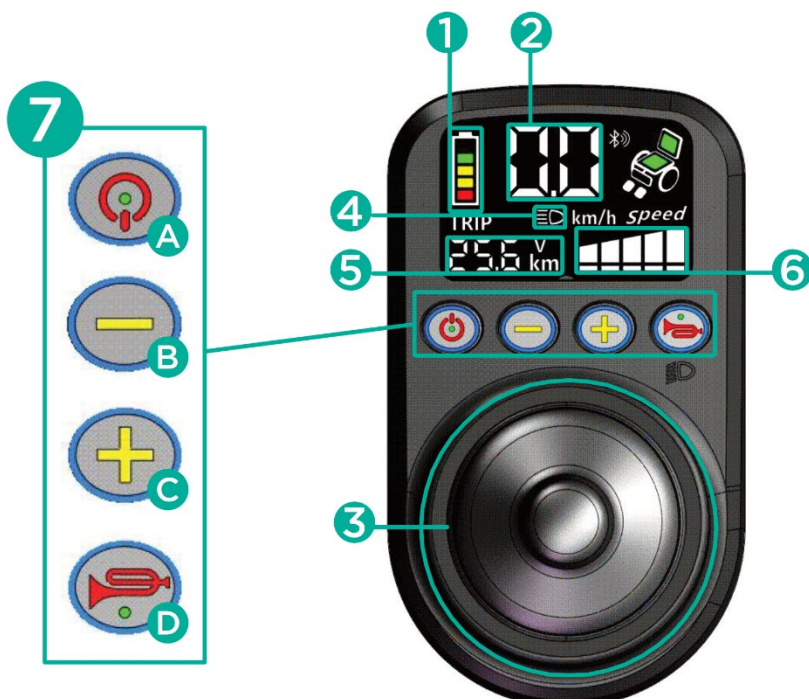
Product code	GO-TIM
Type	for indoor and outdoor use
Overall height	95 cm
Overall width	57 cm
Overall length	96 cm
Seat height from the ground	48 cm (52 cm with cushion)
Backrest height	43 cm
Seat cushion dimensions	41 x 42 x 4 cm
Folded dimensions	72 x 30 x 73 cm
Maximum capacity	136 kg
Battery weight	2,00 kg
Battery	lithium, 24 [V] x 12 [Ah]
Battery charging time	6-8 hours
Wheelchair weight	14,70 kg
Wheels	front: solid, swivel 7" rear: solid 11"
Maximum speed	≤ 6 [km/h]
Braking distance on horizontal surface	≤1,5 [m]
Maximum range	18 [km]
Max. angle of elevation	≤6°
Min. turning radius	1,2 [m]
Safety Classification	type B
Ingress protection level for engine	IP54
Engine power level	2 x 150 [W]

D. Assembly and adjustment

- **Unfolding:** Remove the wheelchair from the box. Place it on the ground. Grab the wheelchair by the armrests and unfold it. Then, press down on the sides of the seat so that all the elements fit together. Lift the backrest by the handles until you hear a “click” sound. Attach the upholstery and lower footrests/anti-tip wheels. Mount control module by sliding it into opening under armrest and securing with adjustment knob (5). Turn on the wheelchair by pressing the power on/off button located on the module.
- **Folding:** To fold the wheelchair, repeat the steps from *Unfolding* in reverse order. **Note!** Before transporting wheelchair, disconnect/turn off battery.
- **Control module adjustment:** Adjustment should be done while seated. Untighten the safety buckle and loosen the adjustment knob (5), adjusting the position or length of the module. Tighten the screw and buckle afterwards.
- **Armrest tilting:** Press the red buttons located on the outside of the wheelchair to lift the armrests. To lower the armrest, pull them down until you hear a “click” sound.
- **Removing/Installing the battery:** To remove the battery, disconnect cable connecting battery to engines (8), press the armrest flip-up button (6) and slowly pull up the battery by the handle. To install battery again, repeat the steps.
- **Anti-tip wheels folding:** Anti-tip wheels can be folded for transportation and storage. Pull lightly a lock (17) and move wheels upward till pin locks in place.

5. OPERATION

A. CONTROL MODULE FUNCTIONS

1. Battery level indicator

2. Current speed display

3. Joystick

4. Flashlight on/off

5. Driving distance

6. Speed levels

7. Operation buttons:

A. Power on/off

B. Speed level down (press)
Warning function (hold)

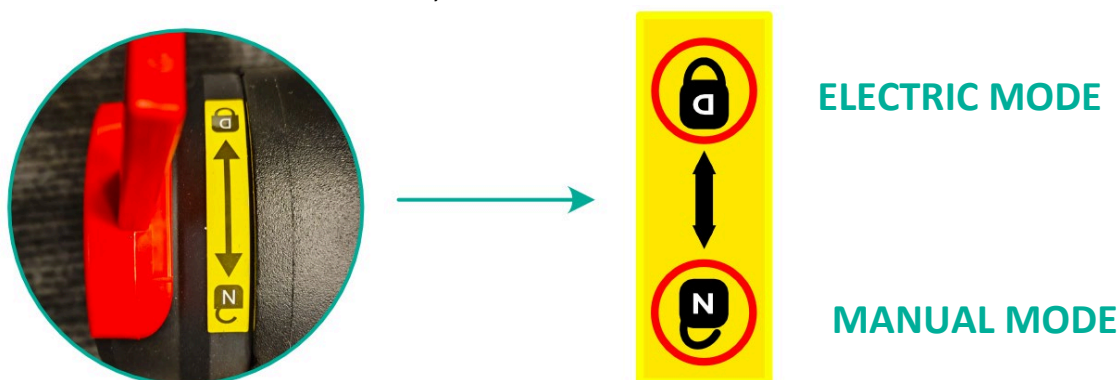
C. Speed level up (press)
Voice assistant on/off
(hold)

D. Flashlight on/off (press)
horn (hold)

Functions	
Joystick	Used to control the speed and direction. The further the joystick is pushed away from the middle, the faster the wheelchair will move. It will lose speed automatically after the joystick is released.
Speed levels indicator	1 – slowest, 5 – fastest.
Voice assistant off/on	Hold button for more than 5 seconds to turn on/off the voice assistant.
Warning function	Hold button to activate warning function. To turn it off, press any other button.
Driving distance	Shows traveled distance when wheelchair is turned on. It resets after turning off the control module.

B. MANUAL AND ELECTRIC MODE SWITCHING

To use the wheelchair in electric mode, move up the left and right engine levers located at the back. To use it in manual mode, move the levers down.

C. CHARGER

- Set includes a charger dedicated to the battery.
- The battery charging time is approx. 6-8 hours.
- Battery can be charged directly or via control module.
- **Important:** when charging the battery through the port under the control module, switch off the module first.
- **Battery LED indicator:**
 - Green LED light** – battery is not connected to the charger or is fully charged.
 - Red LED light** – battery is charging/isn't charged yet

Note! After the charging is finished, the charger turns off to avoid overcharging.

- Check if circuit connection is correct properly and not damaged.
- **Do not** modify charging circuit connection.
- To avoid fire or burnout, avoid unplugging charging circuit while battery is charging.
- To check the current battery level on control module, plug off the charge first.

D. DIRECTIONS FOR USE

- Select the desired driving mode.

Note! Do not switch mode while driving on elevation.

- Press the power button to check the proper functioning of brakes. Brakes are working correctly (electric mode) if you cannot easily push the wheelchair. Otherwise, contact the Seller.
- Tilt the footplates and sit comfortably on the seat. Lower the footplates back to horizontal position.

Note! Do not use the footrest as support when sitting down or standing up from the wheelchair. This may cause the wheelchair to tilt and fall.

- Sit on the wheelchair and turn the control panel on. The LCD panel will light up. Joystick should be in neutral position.
- Take caution when using it for the first time. Control the speed and direction carefully. The electromagnetic brake will get released only after pushing the joystick (electric mode). To speed up, push the joystick further away from its neutral position. To slow down – do the opposite.
- To stop the wheelchair, release the joystick. Pressing the power button will also stop the wheelchair.
- Using dedicated buttons on the control panel, adjust the wheelchair's speed according to user's physical state and road conditions.
- Wheelchair should only be used on flat surfaces. Driving on muddy or uneven surfaces can cause control system damage.

6. MAINTENANCE INFORMATION

1. Wheelchair cleaning and maintenance

- It is recommended to regularly inspect the construction and its elements (wheels, frames, cables, control module, brakes etc.).
- Any used or damaged wheelchair parts need to be replaced (in such case, please contact the Seller).
- Frame should be cleaned with dry cloth to avoid any electrical shorts.

2. Battery maintenance and cleaning:

- 1) Regularly check the battery level displayed. If the level is very low, charge it immediately to avoid a voltage drop (which may affect battery life).
- 2) If wheelchair isn't used for longer period of time, charge battery at least once a month.
- 3) Battery has a certain lifespan. Exchange your battery if it wears out.
- 4) Battery should be kept clean and away from kids. Avoid causing any shock or damage to the battery.
- 5) Charge it regularly to expand its lifespan.
- 6) Used batteries should be utilized according to local environmental regulations.

3. Storage and transport conditions

A. TRANSPORTATION

- Take caution during transportation. Do not throw, turn over or put pressure on the product and its packaging.
- For wheelchair folding, check,, Assembly and adjustment”.

B. STORAGE

Wheelchair should be kept in dry and well-ventilated space. Avoid: high temperature and its rapid changes, acids, alkalis and corrosive chemicals.

C. STORING CONDITIONS

Temperature: -40[°C] ~+55[°C]

Humidity: ≤80%

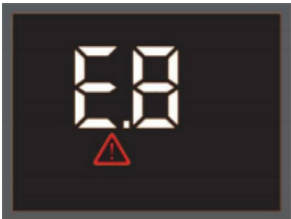
Air pressure: 86 ~ 106 kPA

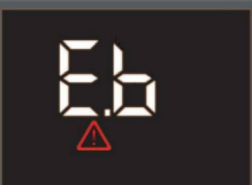
4. Troubleshooting

- In the event of a malfunction, the control module will indicate it with voice notification, sound and warning appearing on the screen. The number displayed on the screen will help to identify the problem.

Issue indicator	Cause	Possible solution
Module doesn't turn ON after pressing power button	Issue with battery or control module	A. Check if battery is charged / power button is on / power cable is correctly mounted / isn't damaged. B. Contact the seller to replace damaged parts.
After turning on, the warning appears on screen, with ring alarm and voice notification. 	Joystick isn't in the center position	Fix the joystick to the center position and reset controller by turning it off. Turn it ON once again – if issue still appears, contact seller for controller calibration.

After turning on, the warning appears on screen, with ring alarm and voice notification. 	The battery level is low.	Charge the battery for minimum 30 minutes before turning it on.
After turning on, the warning appears on screen, with ring alarm and voice notification. 	Left motor is disconnected or damaged	Check motor connection cables. Swap right motor cable with left cable. If the indicator on panel changes to E4 – check if cables are loose or motors are damaged. If E2 code still shows up afterwards, contact seller service.
After turning on, the warning appears on screen, with ring alarm and voice notification 	Left brake is in manual mode / is damaged.	A. Switch left brake to electric mode. B. Check if engines or power cords are damaged. C. Pull off left brake connecting cable and connect to right brake connector. Do the same for right brake connector. If problem still appears, then it's related to the controller. Contact seller for more information.
After turning on, the warning appears on screen, with ring alarm and voice notification 	Right motor is disconnected or damaged	Check motor connection cables. Swap right motor cable with left cable. If the indicator on panel changes to E4 – check if cables are loose or motors are damaged. If E2 code still shows up afterwards, contact seller service.

After turning it on, the warning appears on screen, with ring alarm and voice notification 	Right brake is on manual mode / is damaged.	A. Switch right brake to electric mode. B. Check if engines or power cords are damaged. If yes – contact seller. C. Pull off right brake connecting cable and connect to right brake connector. Do the same for right brake connector. If problem still appears then it's related to the controller. Contact seller for more information.
After turning it on, the warning appears on screen, with ring alarm and voice notification 	Internal controller fault	Contact seller for repair or replacement.
After turning it on, the warning appears on screen, with ring alarm and voice notification 	Error between upper and lower controller (control module and control box)	Contact seller for more information
After turning on, the red warning appears on screen, with ring alarm and voice notification. 	Right/Left Hall sensor issue (brushless motor)	Swap motor cables with each other. If any of the two issues still appear on the panel, contact seller for further information.

		
After turning on, the red warning appears on screen, with ring alarm and voice notification. 	Brake output or short circuit issue	Gently push the joystick. If the wheelchair does not move and you can't hear "click" sound made by electromagnetic brakes, brake or circuit are damaged. Contact seller for repair or replacement.

5. Warranty information

- All products distributed by our company are covered by a warranty, the terms of which are described in the warranty card available on our website. Please contact the retailer from whom you purchased the product. Please, keep in mind that for warranty purposes, the proof of purchase (receipt or invoice) should be retained. As a user-friendly company, we provide only tested products – in terms of material, quality and functionality. If you have any questions regarding service procedures, please contact us.
- **Warranty does not cover:**
 1. **wearable parts:** upholstery, tires, armrests, battery etc.
 2. **damages resulting from:** self-repairs, changes in product construction, dropping, transportation, inappropriate usage, failure to follow instructions, natural disasters (fire, earthquake, flood).

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Zhejiang Feilechi Medical Equipment Co., Ltd.
SUNGO Cert GmbH